

HUMAN SOLUBLE ST2/IL-1R4 ELISA KIT

FOR THE QUANTITATIVE DETERMINATION OF HUMAN ST2 CONCENTRATIONS IN CELL CULTURE SUPERNATES, SERUM AND EDTA PLASMA



ALWAYS REFER TO LOT SPECIFIC PROTOCOL PROVIDED WITH EACH KIT FOR INSTRUCTIONS. PROTOCOL MUST BE READ BEFORE USING THIS PRODUCT.

FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.

PURCHASE INFORMATION:

| ELISA NAME | HUMAN SOLUBLE ST2/IL-1R4 ELISA |
|---|---|
| Catalog No. | SK00120-06 |
| Lot No. | |
| Formulation | 96 T |
| Standard Range | 12.5 - 800 pg/mL |
| Sensitivity | 5 pg/mL |
| Sample Volume | 100 µL |
| Sample Type | Serum, EDTA Plasma, Cell Culture Supernates |
| Dilution Factor | 8 (Optimal dilutions should be determined by each laboratory for each application) |
| Specificity | Human ST2 |
| Calibration | Human ST2 recombinant |
| Intra-assay Precision | 4 - 6% |
| Inter-assay Precision | 8 - 10% |
| Storage | 2 - 8° C |
| This kit contains sufficient materials to run 35 samples duplicated provided that assay is run according to protocol. | |

Order Contact:

AVISCIERA BIOSCIENCE, INC.
2348 Walsh Ave., Suite C
Santa Clara, CA 95051
USA

Tel: (408) 982 0300

Fax: (408) 982 0301

Email: Sales@AvisceraBioscience.com

Info@AvisceraBioscience.com

www.AvisceraBioscience.com

DESCRIPTION

This Human soluble ST2 ELISA Kit contains the necessary components required for the quantitative measurement of recombinant and/or natural human ST2 from cell culture supernates, serum and plasma in a sandwich ELISA format.

This immunoassay contains recombinant human ST2 and antibodies raised against this protein. Results from this immunoassay have shown to accurately quantify recombinant and natural ST2 samples.

ASSAY OVERVIEW

This assay employs the quantitative sandwich ELISA format. The plate is pre-coated with an antibody specific for human ST2. The capture antibody can bind to the human ST2 in the standard and samples. After washing the plate of any unbound substances, an antibody-HRP conjugate against human ST2 is added to the wells. After the last wash to remove any unbound enzyme, a substrate solution is added to the wells and color develops in direct proportion to the amount of human ST2 bound in the standard solutions or samples. A standard curve can be established and sample values can be read off the standard curve.

PROCEDURAL LIMITATIONS

_FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.

_This ELISA kit should not be used beyond the expiration date on the kit label.

_Do not mix reagents with those from other lots or sources.

_It is important that the Dilution Buffer selected for the standard curve be consistent with the samples being assayed.

_Each laboratory must determine the optimal dilution factors for the samples being assayed with a pretest. If samples generate values that are not within the dynamic range of the standard curve, further concentrate or dilute the samples as required with Dilution Buffer and repeat the assay.

_Any modifications in buffers, pipetting technique, washing technique, incubation time or temperature, as well as kit age can cause a change in signal.

_Not all interfering factors have been tested in the immunoassay, therefore the possibility of interference cannot be excluded.

COMPONENTS PROVIDED

| Description | Code | Quantity |
|--|------------------|-----------------|
| ST2 Microplate - 96 well polystyrene microplate (12 strips of 8 wells) coated with an antibody against human ST2. | 120-06-01 | 1 plate |
| ST2 Standard – 3200 pg/vial of recombinant human ST2 in a buffered protein base with preservative; lyophilized. | 120-06-02 | 1 vial |
| Detection Antibody-HRP Conjugate – 105 µL/vial of 100-fold concentrated solution of antibody conjugated to HRP against human ST2. | 120-06-03 | 1 vial |
| Positive Control – one vial of recombinant human ST2; lyophilized. | 120-06-04 | 1 vial |
| Dilution Buffer – 60 mL of buffered protein based solution with preservative. | DB10 | 1 bottle |
| Wash Buffer - 50 mL of 10-fold concentrated buffered surfactant, with preservative. | WB01 | 1 bottle |
| Substrate Solution - 11 mL of TMB substrate solution. | TMB01 | 1 bottle |
| Stop Solution - 11 mL of 0.5M HCl. | S-STOP | 1 bottle |
| Plate Sealer | EAPS | 1 |
| Plastic Pouch | P01 | 1 |

STORAGE

Unopened Kit: Store at 2 – 8° C for up to 6 months. For longer storage, unopened Standard, Positive Control and Detection Antibody-HRP Conjugate 100-fold concentrated solution should be stored at -20° C or -70° C. Substrate Solution can be stored at 2 – 8° C for up to 6 months (**DO NOT FREEZE and PROTECT FROM LIGHT**). All other components can be stored at 2 – 8° C for up to 6 months. Do not use kit past expiration date.

Opened / Reconstituted Reagents: Reconstituted Standard (stock) solution SHOULD BE STORED at -20° C or -70° C for up to one month.

Microplate Wells: Return unused microplate strips to the plastic pouch with the desiccant pack. Microplate may be stored for up to 6 months at 2 – 8 °C.

ADDITIONAL MATERIALS REQUIRED

- Microplate reader capable of absorbance measurement at 450 nm.
- Microplate shaker (250 – 300 rpm).
- Microplate washer or manifold dispenser.
- 100 mL and 500 mL graduated cylinders.
- Multi-channel Pipette, Pipettes and pipette tips.
- Deionized or distilled water.

PRECAUTION

This kit should be handled by those persons who have been trained in and can follow the principles of good laboratory practice. Wear protective clothing such as laboratory overalls, safety glasses and gloves. Care should be taken while handling solutions in this kit to avoid contact with skin or eyes, especially with the stop solution because it contains diluted hydrochloric acid. Wash immediately with water in case of contact on skin or eyes.

SAMPLE COLLECTION AND STORAGE

Cell Culture Supernates - Remove particulates by centrifugation and assay immediately or aliquot and store samples at $\leq -20^{\circ}\text{C}$. Avoid repeated freeze-thaw cycles.

Serum - Use a serum separator tube (SST) and allow samples to clot for 30 minutes before centrifugation for 15 minutes at 1000 x g. Remove serum and assay immediately or aliquot and store samples at $\leq -20^{\circ}\text{C}$. Avoid repeated freeze-thaw cycles.

Plasma - Collect plasma using EDTA, heparin, or citrate as an anticoagulant. Centrifuge for 15 minutes at 1000 x g within 30 minutes of collection. Assay immediately or aliquot and store samples at $\leq -20^{\circ}\text{C}$. Avoid repeated freeze-thaw cycles.

Optional: Use Aprotinin (enzyme inhibitor) for ALL sample collection to prevent sample degradation. 0.5 TIU per ml of sample solution.

SAMPLE PREPARATION

Serum and plasma samples may require dilution. A suggested 8-fold dilution is 40 μL sample + 280 μL Dilution Buffer.

Optimal dilutions should be determined by each laboratory for each application.

Use polypropylene test tubes.


REAGENT PREPARATION

Bring all reagents to room temperature before use.

Wash Buffer - If crystals have formed in the concentrate, warm to room temperature and mix gently until the crystals have completely dissolved. Dilute 50 mL of Wash Buffer Concentrate into deionized or distilled water (450 mL) to prepare 500 mL of 1x Wash Buffer.

ST2 Standard - Reconstitute the ST2 Standard with 1.0 mL of Dilution Buffer. This reconstitution produces a stock solution of 3200 pg/mL. Allow the stock standard to sit for at least 15 minutes with gentle agitation until completely dissolved prior to making standard dilutions (see below). Mix each tube thoroughly before the next transfer. The **800 pg/mL** standard serves as the high standard. The Dilution Buffer serves as the zero standard (0 pg/mL).

| Tube | Standard | Dilution Buffer | Concentration |
|-------|----------------------------|--------------------|---------------|
| Stock | Powder | 1000 μL | 3200 pg/ml |
| # 1 | 125 μL of stock | 375 μL | 800 pg/ml |
| # 2 | 250 μL of 1 | 250 μL | 400 pg/ml |
| # 3 | 250 μL of 2 | 250 μL | 200 pg/ml |
| # 4 | 250 μL of 3 | 250 μL | 100 pg/ml |
| # 5 | 250 μL of 4 | 250 μL | 50 pg/ml |
| # 6 | 250 μL of 5 | 250 μL | 25 pg/ml |
| # 7 | 250 μL of 6 | 250 μL | 12.5 pg/ml |

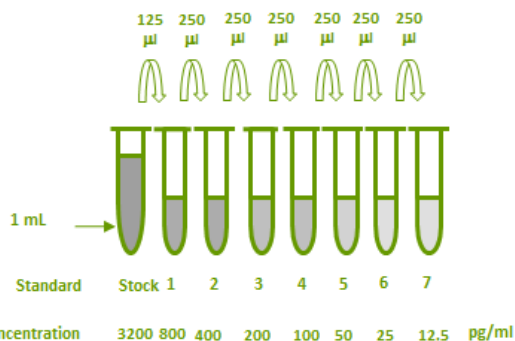

Positive Control - Reconstitute the Positive Control with 1.0 mL of Dilution Buffer. **NOTE:** Positive control could be reused within a few days if stored at $-20^{\circ}\text{C} \sim -70^{\circ}\text{C}$.

Detection Antibody-HRP Conjugate - Pipette 10.395 mL of Dilution Buffer into a 15 mL centrifuge tube and transfer 105 μL of 100-fold concentrated stock solution to prepare working solution. **Note:** 1x working solution of Detection Antibody-HRP conjugate should be used within a few days (**protect from light**). **DO NOT FREEZE.**

ELISA PROTOCOL

Bring all reagents and samples to room temperature before the start of the assay. Blank, standard dilutions, positive control and samples should be assayed in duplicate. ELISA Protocol may need further optimization.

1. Prepare all reagents and working standards as directed in the previous sections.
2. Remove excess microplate strips from the plate frame, return them to the plastic pouch with the desiccant pack.
3. Add 100 μL of Dilution Buffer to Blank wells.
4. Add 100 μL of Standard dilutions, sample, or positive control per well. Cover with plate sealer. Incubate for 2 hours on microplate shaker at room temperature.
5. Aspirate each well and wash, repeating the process three times for a total of four washes. Wash by filling each well with 1x Wash Buffer (300 μL) using a squirt bottle, manifold dispenser, or autowasher. Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Buffer by aspirating or decanting. Invert the plate and blot it against clean paper towels.
6. Add 100 μL of 1x Detection Antibody-HRP conjugate working solution to each well. Cover with plate sealer. Incubate for 1 hour on microplate shaker at room temperature. **Protect from light.**
7. Repeat the aspiration/wash as in step 5.
8. Add 100 μL of Substrate Solution to each well. Incubate for 3-7 minutes on microplate shaker at room temperature. **Protect from light.**
9. Add 100 μL of Stop Solution to each well. The color in the wells should change from blue to yellow. If the color in the wells is green, or if the color change does not appear uniform, gently tap the plate to ensure thorough mixing.
10. Determine the optical density of each well within 15 minutes, using a microplate reader set to 450 nm.

CALCULATION OF RESULTS

Create a standard curve by plotting the log of the known concentrations of the standard dilutions (x-axis) versus the log of its corresponding O.D. (y-axis) and draw the best fit line through the points. It is recommended to use computer software capable of generating a log-log curve fit to more accurately quantify the standard dilutions.

If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.

TYPICAL DATA

This standard curve is provided for demonstration only. A new standard curve should be generated for each set of samples assayed.







| Standard (pg/mL) | Average OD450 (Corrected) |
|------------------|---------------------------|
| Blank | 0 (0.072) |
| 12.5 | 0.024 |
| 25 | 0.047 |
| 50 | 0.112 |
| 100 | 0.216 |
| 200 | 0.436 |
| 400 | 0.851 |
| 800 | 1.641 |

- Lot No.:
- Positive Control:

SPECIFICITY

| PROTEINS | CROSS-REACTIVITY |
|-------------------------|------------------|
| Human ST2 | 100 |
| Human IL-1R1/Fc Chimera | 0 |
| Human IL-1R2/Fc Chimera | 0 |
| Human IL-1RN | 0 |
| Human IL-1 α | 0 |
| Human IL-1 β | 0 |

SUMMARY OF ASSAY PROCEDURE

| |
|--|
| PREPARE REAGENTS, SAMPLES AND STANDARDS |
|  |
| Add 100 µl of standard dilutions, samples, or positive control to each well. Incubate 2 hours on the plate shaker at RT. |
|  |
| Aspirate and wash 4 times. |
|  |
| Add 100 µl per well 1x Detection Antibody-HRP working solution to each well. Incubate 1 hour on the plate shaker at RT. Protect from light. |
|  |
| Aspirate and wash 4 times. |
|  |
| Add 100 µl Substrate Solution to each well. Incubate 3-7 min on the plate shaker at RT. Protect from light. |
|  |
| Add 100 µl Stop Solution to each well. Read 450nm within 15 min. |