

## HUMAN INTERLEUKIN 10 (IL-10) ELISA KIT

FOR THE QUANTITATIVE DETERMINATION OF HUMAN IL-10 CONCENTRATIONS IN SERUM AND PLASMA



ALWAYS REFER TO LOT SPECIFIC PROTOCOL PROVIDED WITH EACH KIT FOR INSTRUCTIONS. PROTOCOL MUST BE READ BEFORE USING THIS PRODUCT.

FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.

### PRODUCT INFORMATION:

ELISA NAME	HUMAN INTERLEUKIN 10 ELISA
Catalog No.	SK00745-01
Lot No.	
Formulation	96 T
Standard Range	4.6 - 300 pg/mL
Sensitivity	2 pg/mL
Sample Volume	100 µL
Dilution Factor	<b>Optimal dilutions should be determined by each laboratory for each application</b>
Sample Type	Serum, Plasma
Specificity	Human IL-10
Calibration	Human IL-10 recombinant
Intra-assay Precision	4 - 10%
Inter-assay Precision	8 - 12%
Storage	2 - 8° C
This kit contains sufficient materials to run 35 samples duplicated provided that assay is run according to protocol.	

### Order Contact:

**AVISCERA BIOSCIENCE, INC.**  
 2348 Walsh Ave., Suite C  
 Santa Clara, CA 95051  
 USA

Tel: (408) 982 0300

Fax: (408) 982 0301

Email: [Sales@AvisceraBioscience.com](mailto:Sales@AvisceraBioscience.com)

[Info@AvisceraBioscience.com](mailto:Info@AvisceraBioscience.com)

[www.AvisceraBioscience.com](http://www.AvisceraBioscience.com)

## DESCRIPTION

This Human Interleukin 10 (IL-10) ELISA Kit contains the necessary components required for the quantitative measurement of recombinant and/or natural human IL-10 from serum and plasma in a sandwich ELISA format.

This immunoassay contains recombinant human IL-10 and antibodies raised against this protein. Results from this immunoassay have shown to accurately quantify recombinant and natural IL-10 samples.

## ASSAY OVERVIEW

This assay employs the quantitative sandwich ELISA format. The plate is pre-coated with an antibody specific for human IL-10. The capture antibody can bind to the human IL-10 in the standard and samples. After washing the plate of any unbound substances, a biotinylated antibody against human IL-10 is added to the wells. After another washing of the plate, Avidin-HRP Conjugate is added. After the last wash to remove any unbound enzyme, a substrate solution (TMB) is added to the wells and color develops in direct proportion to the amount of human IL-10 bound in the standard solutions or samples. A standard curve can be established and sample values can be read off the standard curve.

## PROCEDURAL LIMITATIONS

\_FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.

\_This ELISA kit should not be used beyond the expiration date on the kit label.

\_Do not mix reagents with those from other lots or sources.

\_It is important that the Dilution Buffer selected for the standard curve be consistent with the samples being assayed.

\_Each laboratory must determine the optimal dilution factors for the samples being assayed with a pretest. If samples generate values that are not within the dynamic range of the standard curve, further concentrate or dilute the samples as required with Dilution Buffer and repeat the assay.

\_Any modifications in buffers, pipetting technique, washing technique, incubation time or temperature, as well as kit age can cause a change in signal.

\_Not all interfering factors have been tested in the immunoassay, therefore the possibility of interference cannot be excluded.

## COMPONENTS PROVIDED

DESCRIPTION	CODE	QUANTITY
<b>IL-10 Microplate</b> - 96 well polystyrene microplate coated with an antibody against IL-10.	<b>745-01-01</b>	<b>1 plate</b>
<b>IL-10 Standard</b> – 750 pg/vial of recombinant human IL-10 in a buffered protein base with preservative; lyophilized.	<b>745-01-02</b>	<b>1 vial</b>
<b>Detection Antibody Concentrate</b> – 1.05 mL/vial, 10-fold concentrate of biotinylated antibody against IL-10 with preservative; lyophilized.	<b>745-01-03</b>	<b>1 vial</b>
<b>Positive Control</b> - one vial of recombinant human IL-10; lyophilized.	<b>745-01-04</b>	<b>1 vial</b>
<b>Avidin-HRP Conjugate</b> - 50 µL/vial, 250-fold concentrated solution of Avidin conjugate to HRP.	<b>AVHRP</b>	<b>1 vial</b>
<b>Dilution Buffer</b> – 60 mL of buffered protein based solution with preservative.	<b>DB07</b>	<b>1 bottle</b>
<b>Wash Buffer</b> - 50 mL of 10-fold concentrated buffered surfactant with preservative.	<b>WB01</b>	<b>1 bottle</b>
<b>TMB Substrate Solution</b> - 11 mL of TMB substrate solution.	<b>TMB01</b>	<b>1 bottle</b>
<b>Stop Solution</b> – 11 mL of 0.5M HCl	<b>S-STOP</b>	<b>1 bottle</b>
<b>Plate Sealer</b>	<b>EAPS</b>	<b>1 piece</b>
<b>Plastic Pouch</b>	<b>P01</b>	<b>1 piece</b>

## STORAGE

**Unopened Kit:** Store at 2 - 8° C for up to 8 months. For longer storage, unopened Standard, Positive Control and Detection Antibody Concentrate should be stored at -20° C or -70° C. Do not use kit past expiration date.

**Opened / Reconstituted Reagents:** Reconstituted Standard (stock) solution and Detection Antibody concentrated solution SHOULD BE STORED at -20° C or -70° C for up to one month. Avidin-HRP Conjugate 250-fold concentrated solution and other components may be stored at 2 - 8° C for up to 8 months.

**Microplate Wells:** Return unused wells to the plastic pouch with the desiccant pack. Microplate may be stored for up to 6 months at 2 - 8° C after opening.

**ADDITIONAL MATERIALS REQUIRED**

- Microplate reader capable of absorbance measurement at 450 nm.
- Microplate shaker (250 – 300 rpm).
- Microplate washer or manifold dispenser.
- 100 mL and 500 mL graduated cylinders.
- Multi-channel Pipette, Pipettes and pipette tips.
- Deionized or distilled water.

**SAMPLE COLLECTION AND STORAGE**

**Serum** - Use a serum separator tube (SST) and allow samples to clot for 30 minutes before centrifugation for 15 minutes at 1000 x g. Remove serum and assay immediately or aliquot and store samples at ≤ -20° C. Avoid repeated freeze-thaw cycles.

**Plasma** - Collect plasma using EDTA, heparin, or citrate as an anticoagulant. Centrifuge for 15 minutes at 1000 x g within 30 minutes of collection. Assay immediately or aliquot and store samples at ≤ -20° C. Avoid repeated freeze-thaw cycles.

**SAMPLE PREPARATION**

Serum and Plasma samples do not require dilutions.

**Optimal dilutions should be determined by each laboratory for each application with a pretest. Use polypropylene test tubes.**

**REAGENT PREPARATION**


**Bring all reagents to room temperature before use.**

**Wash Buffer** - If crystals have formed in the concentrate, warm to room temperature and mix gently until the crystals have completely dissolved. Dilute 50 mL of Wash Buffer Concentrate into deionized or distilled water (450 mL) to prepare 500 mL of 1x Wash Buffer.

**IL-10 Standard** - Reconstitute the IL-10 standard with 1.0 mL of Dilution Buffer. This reconstitution produces a stock solution of 750 pg/mL. Allow the standard to sit for a minimum of 15 minutes with gentle agitation prior to making dilutions. Pipette 300 µL of Dilution Buffer into tube #1 and transfer 200 µL of “750 pg/mL” standard solution to make the high standard 300 pg/mL. Pipette 250 µL of Dilution Buffer into tubes #2 to #7. Use the “300

pg/mL” standard solution to produce a dilution series (below). Mix each tube thoroughly before the next transfer. The **300 pg/mL** standard serves as the high standard. The Dilution Buffer serves as the zero standard (0 pg/mL).

TUBE	STANDARD	DILUTION BUFFER	CONCENTRATION
Stock	Powder	1000 µl	750 pg/ml
# 1	200 µl of stock	300 µl	300 pg/ml
# 2	250 µl of 1	250 µl	150 pg/ml
# 3	250 µl of 2	250 µl	75 pg/ml
# 4	250 µl of 3	250 µl	37.5 pg/ml
# 5	250 µl of 4	250 µl	18.75 pg/ml
# 6	250 µl of 5	250 µl	9.375 pg/ml
# 7	250 µl of 6	250 µl	4.688 pg/ml

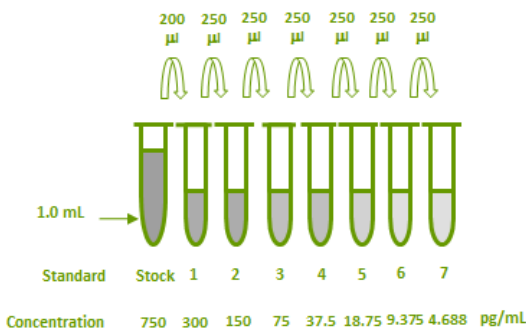

**Positive Control** - Reconstitute the Positive Control with 1.0 mL of Dilution Buffer. Positive Control could be reused within a few days if stored at -20° C or -70° C.

**Detection Antibody Concentrate** - Reconstitute the Detection Antibody Concentrate with 1.05 mL of Dilution Buffer to produce a 10-fold concentrated stock solution. Pipette 9.45 mL of Dilution Buffer into a 15 mL centrifuge tube and transfer 1.05 mL of 10-fold concentrated stock solution to prepare working solution.

**Avidin-HRP Conjugate** - Pipette 11.952 mL of Dilution Buffer into a 15 mL centrifuge tube and transfer 48 µL of 250-fold concentrated stock solution to prepare working solution. **Note:** 1x working solution of Avidin-HRP should be used within a few days (**protect from light**). **DO NOT FREEZE.**

## ELISA PROTOCOL

Bring all reagents and samples to room temperature before the start of the assay. Blank, standard dilutions, positive control and samples should be assayed in duplicate. ELISA Protocol may need further optimization.

1. Prepare all reagents and working standards as directed in the previous sections.
2. Remove excess microplate strips from the plate frame, return them to the plastic pouch with the desiccant pack.
3. Add 100  $\mu$ L of **Dilution Buffer** to Blank wells.
4. Add 100  $\mu$ L of **Standard dilutions (#7 - #1)** in reverse order of serial dilution, **sample**, or **positive control** per well. Cover with plate sealer. Incubate for 2 hours on microplate shaker at room temperature.
5. Aspirate each well and wash, repeating the process three times for a total of four washes. Wash by filling each well with **1x Wash Buffer** (300  $\mu$ L) using a squirt bottle, manifold dispenser, or autowasher. Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Buffer by aspirating or decanting. Invert the plate and blot it against clean paper towels.
6. Add 100  $\mu$ L of **Detection Antibody working solution** to each well. Cover with plate sealer. Incubate for 2 hours on microplate shaker at room temperature.
7. Repeat the aspiration/wash as in step 5.
8. Add 100  $\mu$ L of **Avidin-HRP Conjugate working solution** to each well. Incubate for 30 minutes on microplate shaker at room temperature. **Protect from light.**
9. Repeat the aspiration/wash as in step 5.
10. Add 100  $\mu$ L of **Substrate Solution** to each well. Incubate for 8-15 minutes on microplate shaker at room temperature. **Protect from light.**
11. Add 100  $\mu$ L of **Stop Solution** to each well. The color in the wells should change from blue to yellow. If the color in the wells is green, or if the color change does not appear uniform, gently tap the plate to ensure thorough mixing.
12. Determine the optical density of each well within 15 minutes using a microplate reader set to 450 nm.

## CALCULATION OF RESULTS

Create a standard curve by plotting the log of the known concentrations of the standard dilutions (x-axis) versus the log of its corresponding O.D. (y-axis)

and draw the best fit line through the points. It is recommended to use computer software capable of generating a log-log curve fit to more accurately quantify the standard dilutions.

If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.

Calculation of samples with a concentration exceeding that of the high standard may result in inaccurate, low human IL-10 levels. Such samples require further external pre-dilution according to expected human IL-10 values with Dilution Buffer in order to precisely quantify the actual human IL-10 level.

## TYPICAL DATA

This standard curve is provided for demonstration only. A standard curve should be generated for each set of samples assayed. The blank is used as a reference for your blank. The O.D. values for the standards have already been corrected (i.e. minus the blank).

STANDARD (PG/ML)	AVERAGE OD450 (CORRECTED)
Blank	0 (0.098)
4.688	0.052
9.375	0.095
18.75	0.112
37.5	0.226
75	0.373
150	0.714
300	1.660

- **Lot No.:**
- **Positive Control:**

## EXPECTED VALUES









There were no detectable human IL-10 levels found in our samples.

Elevated human IL-10 levels depend on the type of immunological disorder.

## SPECIFICITY

PROTEINS	CROSS-REACTIVITY (%)
Human IL-10	100
Human IL-6	0
Human IL-1 $\beta$	0
Human IL-33	0

**SUMMARY OF ASSAY PROCEDURE**

<b>PREPARE REAGENTS, SAMPLES AND STANDARDS</b>

Add 100 µl of standard dilutions, samples, or positive control to each well. Incubate 2 hours on the plate shaker at RT.

Aspirate and wash 4 times.

Add 100 µl Detection Antibody working solution to each well. Incubate 2 hours on the plate shaker at RT.

Aspirate and wash 4 times.

Add 100 µl Avidin-HRP Conjugate working solution to each well. Incubate 30 minutes on the plate shaker at RT. <b>Protect from light.</b>

Aspirate and wash 4 times.

Add 100 µl Substrate Solution to each well. Incubate 8-15 minutes on the plate shaker at RT. <b>Protect from light.</b>

Add 100 µl Stop Solution to each well. Read at 450 nm within 15 min.